



Uricult[®] Vet

Model Chart



Uricult[®] Vet

Accurate and easy to use bacterial culture method for the diagnosis of urinary tract infections (UTI)

The Uricult Vet dipslide carries two culture media, CLED (green) and MacConkey (reddish-brown) medium for detection of microbes causing urinary tract infections in animals. The CLED is for detection of total bacterial count and MacConkey is for detection of gram-negative microbes.

EASY TO USE Dip – Incubate – Interpret

ON-SITE TEST No need to transport the urine sample

SAFE TRANSPORTATION Maintains integrity of the sample


LONG SHELF-LIFE Can be stored at room temperature for nine months

Make your practice stand out

Uricult Vet offers veterinary clinics a fast and convenient on-site method to identify UTI, with accurate results available in 16–24 hours. This practise helps vets to avoid sending negative samples to reference laboratory.

Further information: www.oriondiagnostica.com

IT'S SO EASY TO USE!

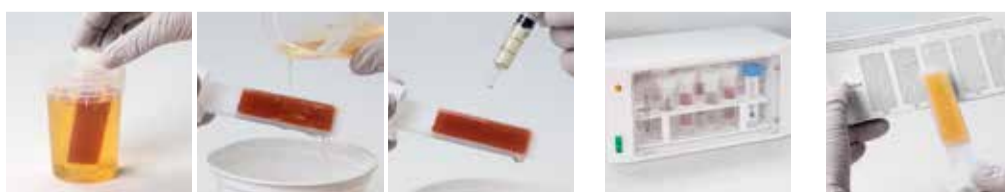

Dip, pour or drop

Incubate

Interpret

For more detailed information, please consult the instructions for use.



Product list	Cat. No.
Uricult Vet 10 tests	06064
Cultura M 230 V	129997 / 05866
Cultura M 110/115 V	132754 / 06266
Cultura M 230 V (with a UK plug)	133192
Cultura M 12 V (with a car plug)	132697
Thermometer, accessory	129981 / 05724
Multirack, accessory	139248

Orion Diagnostica develops, manufactures and markets diagnostic test systems for healthcare professionals with over 40 years of experience. Building well-being is the core of our mission.

Uricult[®] Vet is a registered trademark of Orion Diagnostica Oy.



Orion Diagnostica Oy, Koivu-Mankkaan tie 6 B, P.O. Box 83, FI-02101 Espoo, Finland
Tel. +358 10 4261, Fax +358 10 426 2794, www.oriondiagnostica.com

8079-01EN, 07/2017. International version.

Organisms on media



Colony densities *Escherichia coli* CFU/ml

